

# ET225 Electrophoresis Platform

## User Manual

Version 1.03 March 2012

T. +61 2 9807 8855

F. +61 2 9807 8844

E. [info@edaq.com](mailto:info@edaq.com)

[www.edaq.com](http://www.edaq.com)

© 2010 eDAQ Pty Ltd.

All rights reserved.

---

## Contents

Introduction .....	3
System Description .....	3
Exploded View.....	4
Assembly .....	5
HV Connections.....	7
Micronit Chip Specification.....	8

# ET225 C4D MicroChip Electrophoresis Platform

---

## Introduction

The ET255 platform was designed to allow standard 45mm and 90mm Microchips manufactured by: Micronit Microfluidics (see information at the end of this document) to be used as a complete microchip based Electrophoresis system.

The complete ER255 (and ER455) system is composed of the following elements:

ER225 C4D System that provides the excitation and signal acquisition components of the system.

ER230 HV Sequencer that provides the HV necessary for the Electrophoresis function

ET225 the mechanical platform (adapter) that provides the required facilities to load samples and acquire the data generated.

ET145-4 Four 45mm borosilicate glass TYPE T35100C4D microchips

EC 20 Standard test solutions.

*Disclaimer: This system is supplied as a research grade instrument and eDAQ does not claim that this system is suitable for any particular Electrophoresis application. The user takes full responsibility for the design of experiments and the interpretation of the resulting data.*

## System Description

A fully wired ET225 system is shown below, it consists of:

Five (5) Multicoloured HVcables

Interface cable


Safety Interlock cable

Ground pin



An exploded view of the system is shown below:

## Exploded View

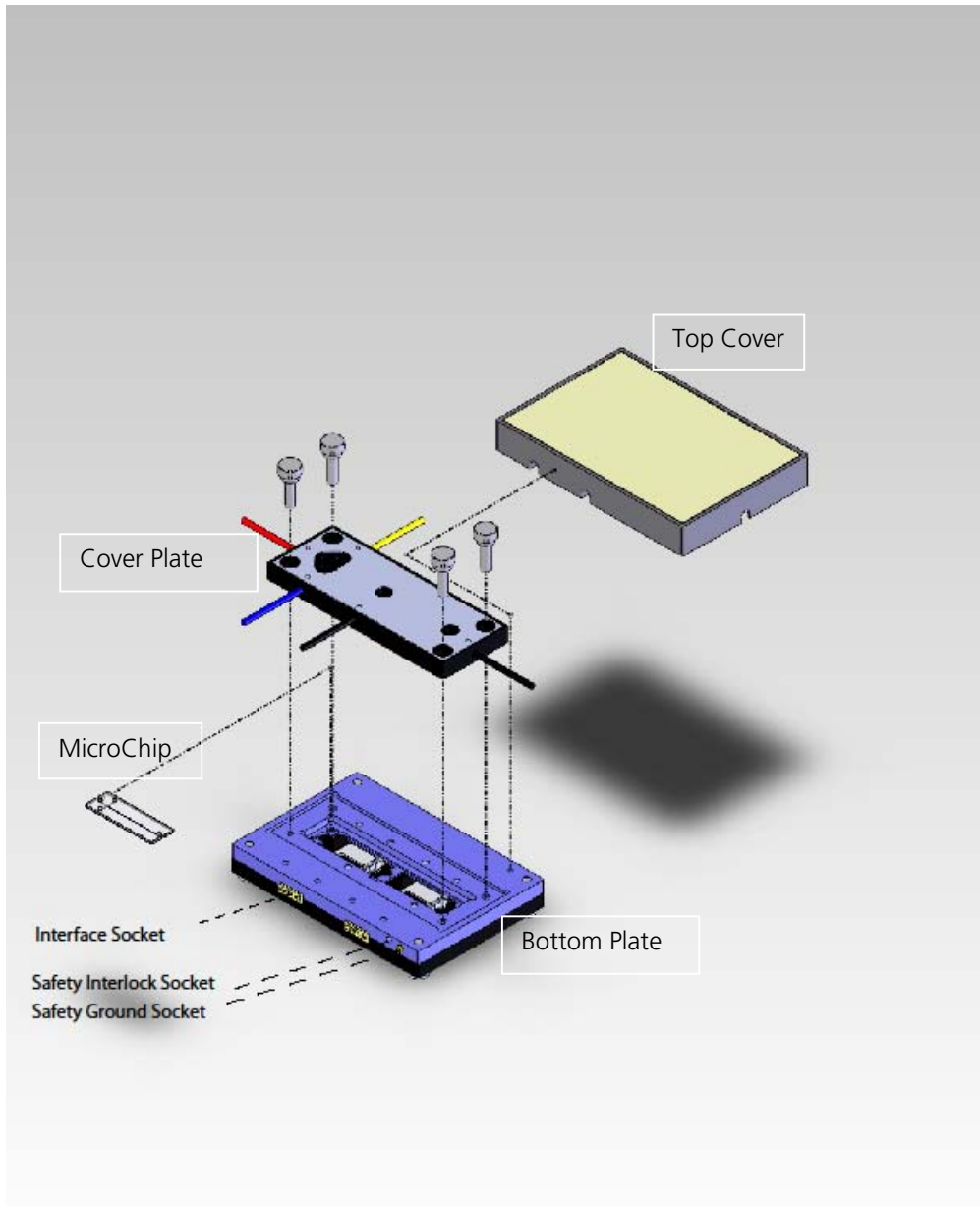

## Assembly

**CAUTION:** Microchips are constructed from glass and should therefore be handled very carefully to prevent breakages. Warranty does not apply to broken or used microchips. Microchips are delivered in a shock proof container.


**Step 1** Connect Interface cable and Safety interlock cables.

**Step 2** Fit ET145 Microchip as shown – ensure that the chip is mounted correctly.



This interlock "Key" has a small flat which should be facing downward when inserting the key in the interlock position. The interlock Key should be fully inserted otherwise it will not operate correctly.

**Step 3** Check that the O-rings are in place before fitting the cover plate – spare O-rings are provided

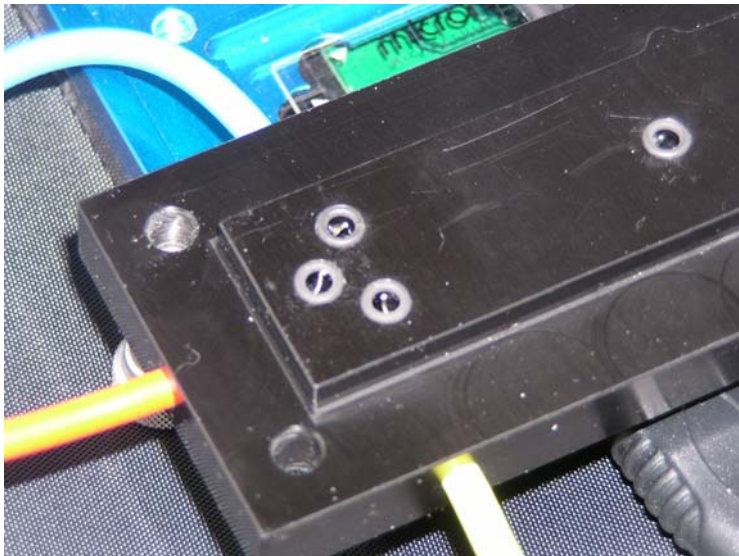


**Step 4** Carefully fit the cover plate using the knurled bolts.

These bolts should be tightened evenly and in turn and **should not be over tightened since this WILL damage the microchips.**

Light finger pressure on the bolts is sufficient to create a good seal between the O-rings and the chip.



The diagram below (NOT TO SCALE) shows diagrammatically the fluid circuit from the Filling well to the microchip fluid channel. The O-ring represents the only contact between the microchip and the cover plate.

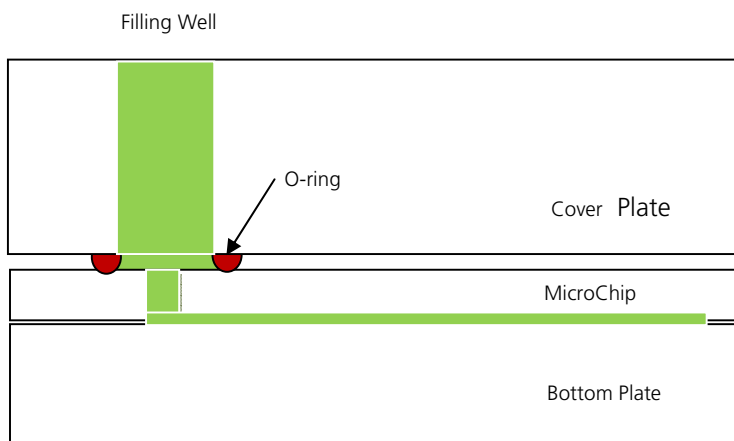

The top cover should now be located, a weak magnetic attraction is provided to keep the Top cover in place and a magnet located in the Top cover over the safety interlock causes the safety switch to be closed thus enabling the safe application of high voltages to the system.


## HV Connections

The Microchip wells are numbered 1 to 4 and since the system is symmetrical different arrangements will provide similar results however it is good practice to set a standard connection configuration and stick to it in order to ensure consistency.

Typically the following HV connections will be made as shown in the diagram below

SAMPLE PORT (1): This will have a high voltage applied during the injection or loading phase and then a HiZ mode during the Separation phase.

**CONNECT BLUE cable to Output 1**

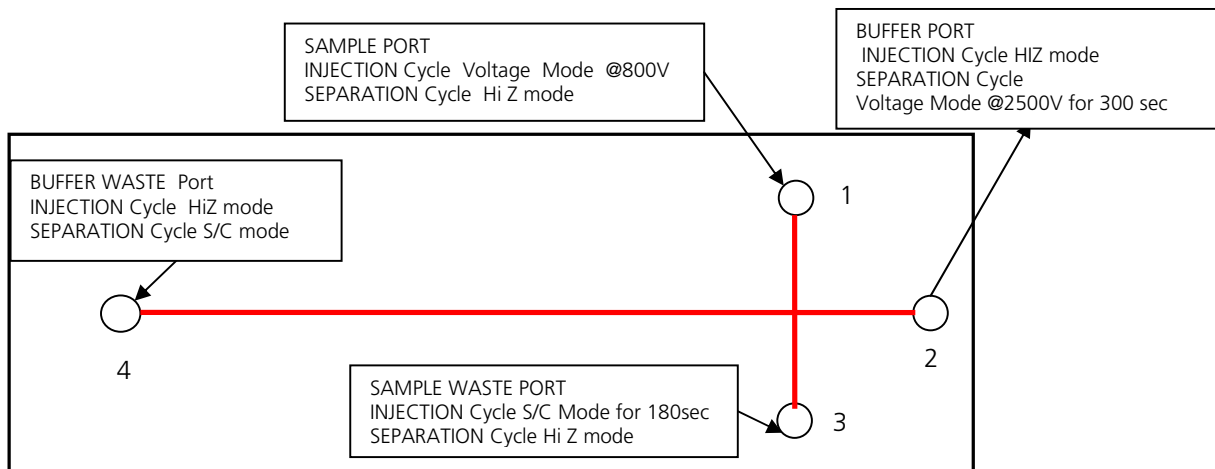

BUFFER PORT (2) Hi Z during injection and a Hi Voltage during separation phase.

**Connect RED cable to Output 2**

SAMPLE WASTE PORT (3) S/C during injection and HiZ during separation phase.

**Connect YELLOW cable to Meter 1**

BUFFER WASTE PORT (4) HiZ during Injection and S/C during separation phase.

**Connect BLACK cable to Meter 2**

Note 1: Two BLACK Cables are provided – Number 4 is used with 45mm microchips and Number 5 is used with 90mm microchips.

Note 2: the detection electronics can be fitted either in the 45mm position or the 90mm position. If both types of microchips will be used routinely then duplicate the detection electronics

## Micronit chip cleaning procedure

INSIDE LIFE SCIENCES

**micronit**  
MICROFLUIDICS

### **Microchip cleaning procedures**

One simple but very effective way to clean a microchip is to flush an alkaline solution through the channels. A solution of 1 M sodium hydroxide in water works well but a lower concentration might also be sufficient. If traces of the cleaning solution remaining inside the chip after cleaning and rinsing with water pose a problem then ammonia can be used instead.

Note that these solutions are caustic and can cause damage to e.g. the polyimide coating of fused silica capillaries. Also the plastic parts of the EOF kit should not be exposed to very alkaline solutions.

In order to aid in the removal of particulate matter a water bath with ultrasonic agitation can be used, preferably while flushing a watery solution through the channels using an Fluidic connect kit.

Glass microchips can be heated (e.g.  $>400^{\circ}\text{C}$ ) causing any organic material adsorbed on the glass surface to degrade.

Concentrated sulphuric acid works well to dissolve organic material such as fibres which are difficult to remove with alkaline solutions, but because of the extremely corrosive nature of the material a cleaning procedure is not so easily implemented.

Visiting address Colosseum 15 • 7521 FV Enschede Postal address P.O. Box 214 • 7500 AE Enschede • The Netherlands  
Phone +31 53 850 6 850 Fax +31 53 850 6 851 E-mail: info@micronit.com Internet: www.micronit.com



## Micronit Chip Specification

See next Page

## Specification document T35100C4D chip

### Chip:

Material: D263 T borosilicate glass

Chip size (w x h): 15.0 x 45.0 +/- 0.3 mm

Thickness: 1.1 +/- 0.050 mm cover and 0.100 +/- 0.015 mm bottom

On-chip reservoir volume: <1  $\mu$ L

### Channels:

Width: 100 +/- 5  $\mu$ m


Depth: 10.0 +/- 0.5  $\mu$ m

Length of the double-T: 100  $\mu$ m between center lines

Channel length (2-4): 40mm

Separation length: 33mm

Cross Arm length: 9mm

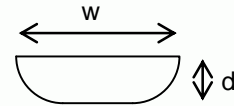

### Electrodes:

Material: platinum, tantalum adhesion layer deposited on the bottom surface

Electrode thickness: 200 nm

Size (wxh): 200 x 500  $\mu$ m

Distance between excitation and detection electrode: 250  $\mu$ m

Distance between electrode pairs: 2.9 mm

